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INSTRUCTION MANUAL EN

ELECTRONIC BABY SCALE
METRICBABY

BEFORE USING THIS DEVICE, PLEASE READ THIS USER MANUAL! FOLLOW ALL SAFETY INSTRUCTIONS TO AVOID DAMAGE CAUSED BY IMPROPER USE! KEEP THIS USER MANUAL FOR FUTURE REFERENCE.


This electronic scale is a device for checking the baby's weight at home, it is especially useful for weighing infants and premature babies, where the weight is an exceptionally important development factor. We do our best efforts to ensure the product with top quality, comfortable and safe during weighing.

THIS DEVICE IS ONLY HOUSEHOLD USE, NO MEDICAL OR COMMERCIAL USE.

SAFETY INSTRUCTIONS

- Ensure to place the scale on flat, dry and hard surface when weighing.
- Use this device for its intended purpose at home.
- Do not load object that over the scale's max capacity: 30kg.
- Avoid shaking or hitting the scale.
- Do not submerge it into water.
- When in storage, do not place anything on the scale.
- The optimal temperature for operating and storing the scale is from 5 to 30 degree celcius.

PRODUCT DESCRIPTION

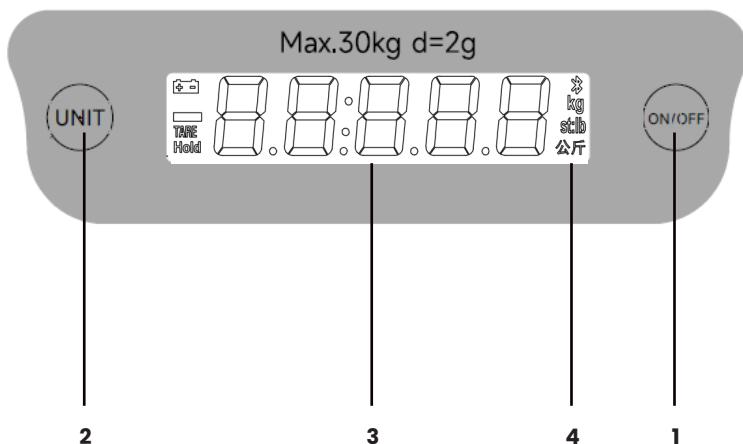1. On/off butoon / tare button	•Press the button to switch on or off the device •While use a towel for more comfortable weighing ,first put the towel on the platform and then press the button to set the display to“0.00”
2. Unit button	•Press the button to convert the units of measurement : gram(g),kilogram(kg),pound(lb)or stone(st:lb)
3. LCD screen	•Display weight
4. G/KG/LB/LB:ST indication	•The unit digits comes on when gram(g), kilogram (kg), pound (lb) or stone (lb:st) is set as unit of measurement.

USER INSTRUCTIONS

Weight

- 1) Please use the scale one the hard & flat surface. Do not weight on carpet or any soft surface as it may affect accuracy.
- 2) Press “ON/OFF” button, the display shows “Hi” ;
- 3) After display “0.00kg”, then ready for weighing ;
- 4) Put baby or other stuff on the scale platform, and LCD display will shows the weight;
- 5) Press “ON/OFF” key to turn off the scale , or scale will be auto off after 30 seconds while in idle.

Auto lock function

When scale has the weight of the baby, do not move the baby off the platform, the result lock on display automatically within 1-2 seconds.

Tare function

It could weight continuously:

- First, put the first baby on the scale platform to get weight, take away the baby,the ba-by's weight still on the LCD.
- Second, press “Tare” button, LCD will turn to zero, then put the next one on the scale for weighing. (this function is repeatable)

Unit conversion

Press “unit” button to convert between“g” “kg” “lb” and “st:lb”.

Overload Indication

When the weight is over the “max capacity” LCD display will show “ 0-Ld”

Low Power Indication

When LCD display shows “LO”, Please replace the new battery with correct polarity.

Battery Installation

- 1) Before using the scale, check the battery compartment and remove any plastic insulation strips
- 2) Insert the batteries (3 x 1.5V AAA) into the battery compartment
- 3) If you are replacing the batteries,make sure the positive and negative contacts are properly aligned.

CLEANING AND STORAGE

- Clean with a soft, slightly damp cloth.
- Use a mild cleaning agent if necessary.
- Do not submerge the scale in water or rinse it under running water.
- Do not use strong detergents, solvents, or corrosive substances.
- All parts must be completely dry before reuse.
- Store in a dry place, away from heat sources.
- Remove the batteries if the scale will not be used for an extended period.

TECHNICAL DATA

CAPACITY	40g ~ 30KG
PRESS ON AND AUTO OFF	yes
DIVISION	2G
4 HIGH PRECISION STRAIN GAUGE SENSORS	yes
TARE FUNCTION	yes
LOW BATTERY AND OVER LOAD INDICATE	yes
AUTO LOCK/HOLD FUNCTION (EASIER WEIGHING OF MOVING BABIES, THE RESULT LOCK ON DISPLAY AUTOMATICALLY WITHIN 1-2 SECONDS)	yes
LARGE LCD	81x16.9mm with backlight
BATTERY	3 x 1.5V AAA batteries
DEVICE DIMENSIONS	700x319x42.5mm

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

DECLARATION OF CONFORMITY

Malis B. Machoński Sp. K. declares that device described in this manual meets all the requirements of the:

(EMC) Directive 2014/30/EU
(RoHS) 2011/65/EU + (EU) 2015/863
(REACH) Regulation (EC) No 1907/2006

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